

## THOMASCOLOR Wax Carver Instruments



### DESCRIPTION

Instruments for wax and composite modelling in esthetic dentistry. The slim, light aluminium handle allows easy handling, and the working tips, made of high-quality stainless steel suitable for use in the oral cavity, guarantee greater durability and minimum wear.

### SPATULAS

#### Thomascolor 1 (yellow) TC1

Sculptor for wax, ceramic and acrylic, with straight handle. Autoclavable. Tips finished in different fine diameters to facilitate waxing.

#### Thomascolor 2 (green) TC2

Sculptor for wax, ceramic and acrylic, with straight handle. Autoclavable. Tips finished in different thicker diameters to facilitate waxing.

#### Thomascolor 3 (blue) TC3


Modeller for wax, ceramic, composite and acrylic. Autoclavable. Special working tips for obturation.

#### Thomascolor 4 (red) TC4

Instrument for modelling wax, ceramic, composite and acrylic. Hardened stainless steel working tips for strength and durability. Two ends with opposite inclination, both with a sharp edge for cutting excess material.

#### Thomascolor 5 (grey) TC5

Instrument for modelling wax, ceramic, composite and acrylic. Two curved ends with spherical tips of 0.45 and 0.60 mm diameter. Hardened stainless steel working tips for greater resistance and duration, specially designed for cutting excess material and for the final profiling of the piece.

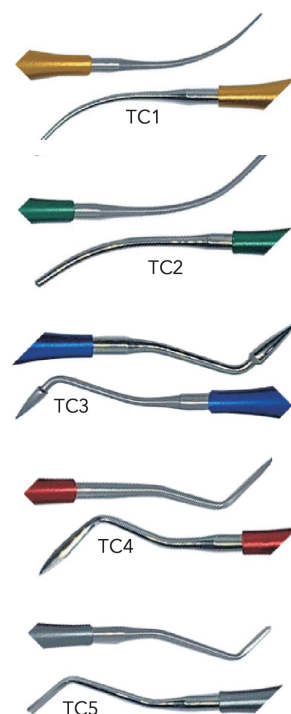

### PRODUCT INFORMATION

#### INTENDED USE

Hand instruments for modelling wax, composite, ceramic and acrylic materials in esthetic dentistry and in the dental laboratory.

#### WARNINGS AND PRECAUTIONS

Check the instrument before each use. Do not use it if a working tip is damaged, bent or loose. The sharp edges of TC4 and TC5 can cut: handle with care.

#### CLEANING AND STERILISATION

Reusable instrument. Clean and sterilise before first use and after each use according to the SAFIDENT Manual Cleaning: steam sterilisation, maximum 150 °C (300 °F).

**Do not place in ultrasonic cleaner!**

#### MATERIAL

Handle: anodised aluminium. Working tips: stainless steel (hardened in TC4 and TC5).

#### MANUFACTURER

SWISS-DS-SAFIDENT, S.L.U. · C/ Ponent 10 · 43700 El Vendrell (Tarragona) · Spain · [info@safident.eu](mailto:info@safident.eu) · [www.safident.eu](http://www.safident.eu)

Class I medical device. For professional dental use only.

